


Midas® III-Plus Automated Stainer

Automated slide stainer for hematology and bacteriology

The Midas® III-Plus Automated Stainer (P/N 64000H) gives you the ability to automate all hematology and bacteriology staining applications. Since the unit is fully programmable, you can duplicate your existing staining techniques and ensure consistent results for your peripheral blood and bone marrow smears as well as your bacteriological smears. MilliporeSigma backs up the Midas® III-Plus Automated Stainer with a one year warranty and some of the best technical and customer service in the business. Don't forget we have a full line of stains (including the Harleco® product line) perfectly suited for use with the Midas® III-Plus Automated Stainer.

Key Features:

- Ideally suited for Wright (P/N 740) & Wright Giemsa (P/N 742) staining of blood and bone marrow smears
- Ease of use for Gram and TB staining applications
- Attractive blue/light blue color for the instrument
- Optional clear colorless enclosure with door

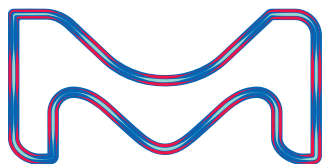

The Cornerstone of Your Laboratory

The Midas® III-Plus Automated Stainer is a fully automated, benchtop stainer with a small footprint. It is designed to perform routine hematology and bacteriology staining effortlessly. It utilizes the dip slide staining technique which ensures uniform staining for every sample. This reduces stain consumption up to 50%.

This versatile automated stainer combines the reliable electronics and mechanics of the original Midas® III Stainer, now with a new optional plastic hood enclosure to make your laboratory a safer environment. Ideally suited for staining up to 300 slides per day, this stainer will become the new workhorse of your hematology/bacteriology lab.

Technical Specifications

Specification	Description
Staining Control	Microprocessor
Available programs	9 user-defined programs (up to 29 steps each), 1 test program
Application time	1 second to 99 minutes
Available stations	6 stations: 4 reagent, 1 drying, 1 rinse (Option to turn drying and rinse stations to reagent)
Stain Vessel	8 vessels & lids
Capacity	300 ml
Slide Carrier	1 20 position slide carrier
Capacity	20 or 30 slides
Rinse Station Flow	Program maintains 500 ml/min increments allowing for precise rinsing
Length	21" (53 cm)
Width	18" (45 cm)
Height	15" (38 cm)
Weight	28 lbs. (13 kg)
Shipping Weight	60 lbs.
UL/cUL approved	

Ordering Information

Product	Order No.	Size	Package
Midas® III-Plus Slide Stainer	640000H-94	1EA	Chipboard Box
Midas® Coverslip Carrier	64053-94	1EA	Chipboard Box
Midas® Slide Carrier 20 slides	64052-94	1EA	Chipboard Box
Midas® Slide Carrier 30 slides	64052E-94	1EA	Chipboard Box
Midas® Vessels and Lids White 8/pk	64054W-95	1EA	Chipboard Box
Midas® Water Bath	64000WB-94	1EA	Chipboard Box
Midas® Water Bath Filters 2/pk	64000WBF-95	1EA	Chipboard Box
Midas® Vessel Rack	RD11953	1EA	Chipboard Box
Midas® III-Plus Operator's Manual	5.06115.0001	1EA	
Midas® Plastic hood enclosure	5.06117.0001	1EA	
Midas® 110V power cord (US)	RW918003	1EA	

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

Indication for Use:
For *in vitro* diagnostic use

Intended User:
For professional use only



To place an order or receive technical assistance

In the U.S. and Canada, call toll-free 1-800-645-5476
For other countries across Europe and the world, please visit: EMDMillipore.com/offices
For Technical Service, please visit: EMDMillipore.com/techservice

EMDMillipore.com/MidasIII-PlusStainer

